



Glucose-6-Phosphate Dehydrogenase Microplate Assay Kit User Manual

Catalog # FTA0045

(Version 1.2D)

Detection and Quantification of Glucose-6-Phosphate
Dehydrogenase (G6PDH) Activity in Urine, Serum, Plasma, Tissue
extracts, Cell lysate, Cell culture media and Other biological fluids
Samples.

For research use only. Not for diagnostic or therapeutic procedures.

I. INTRODUCTION.....	2
II. KIT COMPONENTS.....	3
III. MATERIALS REQUIRED BUT NOT PROVIDED.....	3
IV. SAMPLE PREPARATION.....	4
V. ASSAY PROCEDURE.....	5
VI. CALCULATION.....	6
VII. TYPICAL DATA.....	7
VIII. TECHNICAL SUPPORT.....	7
IX. NOTES.....	7

I. INTRODUCTION

Glucose 6 phosphate dehydrogenase (G6PDH) is a cytosolic enzyme in the pentose phosphate pathway, a metabolic pathway that supplies reducing energy to cells (such as erythrocytes) by maintaining the level of the co-enzyme nicotinamide adenine dinucleotide phosphate (NADPH). The NADPH in turn maintains the level of glutathione in these cells that helps protect the red blood cells against oxidative damage. Of greater quantitative importance is the production of NADPH for tissues actively engaged in biosynthesis of fatty acids and/or isoprenoids, such as the liver, mammary glands, adipose tissue, and the adrenal glands.

The reaction velocity is determined by measuring the increase in absorbance at 340 nm resulting from the reduction of NADP.

II. KIT COMPONENTS

Component	Volume	Storage
96-Well Microplate	1 plate	
Assay Buffer	30 ml x 4	4 °C
Substrate Diluent	20 ml x 1	4 °C
Substrate	Powder x 1	-20 °C
Standard	Powder x 1	-20 °C
Positive Control	Powder x 1	-20 °C
Technical Manual	1 Manual	

Note:

Substrate: add 19 ml Substrate Diluent to dissolve before use.

Standard: add 1 ml distilled water to dissolve before use; then add 0.2 ml into 0.8 ml distilled water, the concentration will be 400 $\mu\text{mol/L}$.

Positive Control: add 1 ml distilled water to dissolve before use, then add 0.1 ml into 0.9 ml distilled water, mix.

III. MATERIALS REQUIRED BUT NOT PROVIDED

1. Microplate reader to read absorbance at 340 nm
2. Distilled water
3. Pipettor, multi-channel pipettor
4. Pipette tips
5. Mortar
6. Centrifuge
7. Timer
8. Ice

IV. SAMPLE PREPARATION

1. For cell and bacteria samples

Collect cell or bacteria into centrifuge tube, discard the supernatant after centrifugation, add 1 ml Assay buffer for 5×10^6 cell or bacteria, sonicate (with power 20%, sonication 3s, interval 10s, repeat 30 times); centrifuged at 8000g 4 °C for 10 minutes, take the supernatant into a new centrifuge tube and keep it on ice for detection.

2. For tissue samples

Weigh out 0.1 g tissue, homogenize with 1 ml Assay buffer on ice, centrifuged at 8000g 4 °C for 10 minutes, take the supernatant into a new centrifuge tube and keep it on ice for detection.

3. For serum or plasma samples

Detect directly.

V. ASSAY PROCEDURE

Warm all reagents to room temperature before use.

Add following reagents in the microplate:

Reagent	Sample	Standard	Blank	Positive Control
Standard	--	200 μ l	--	--
Distilled water	--	--	200 μ l	--
Substrate	190 μ l	--	--	190 μ l
Sample	10 μ l	--	--	--
Positive Control	--	--	--	10 μ l
Mix, measured at 340 nm and record the absorbance of 10th second and 130th second.				

Note:

- 1) Perform 2-fold serial dilutions of the top standards to make the standard curve.
- 2) For unknown samples, we recommend doing a pilot experiment & testing several doses to ensure the readings are within the standard curve range. If the enzyme activity is lower, please add more sample into the reaction system; or increase the reaction time; if the enzyme activity is higher, please dilute the sample, or decrease the reaction time.

VI. CALCULATION

Unit Definition: One unit of G6PDH is defined as the enzyme produce 1 nmol NADPH per minute.

1. According to the protein concentration of sample

$$\begin{aligned} \text{G6PDH (U/mg)} &= (C_{\text{Standard}} \times V_{\text{Standard}}) \times (OD_{\text{Sample}(130\text{S})} - OD_{\text{Sample}(10\text{S})}) / (OD_{\text{Standard}} - \\ &\quad OD_{\text{Blank}}) / (V_{\text{Sample}} \times C_{\text{Protein}}) / T \\ &= 4000 \times (OD_{\text{Sample}(130\text{S})} - OD_{\text{Sample}(10\text{S})}) / (OD_{\text{Standard}} - OD_{\text{Blank}}) / C_{\text{Protein}} \end{aligned}$$

2. According to the weight of sample

$$\begin{aligned} \text{G6PDH (U/g)} &= (C_{\text{Standard}} \times V_{\text{Standard}}) \times (OD_{\text{Sample}(130\text{S})} - OD_{\text{Sample}(10\text{S})}) / (OD_{\text{Standard}} - \\ &\quad OD_{\text{Blank}}) / (V_{\text{Sample}} \times W / V_{\text{Assay}}) / T \\ &= 4000 \times (OD_{\text{Sample}(130\text{S})} - OD_{\text{Sample}(10\text{S})}) / (OD_{\text{Standard}} - OD_{\text{Blank}}) / W \end{aligned}$$

3. According to the quantity of cells or bacteria

$$\begin{aligned} \text{G6PDH (U/10}^4\text{)} &= (C_{\text{Standard}} \times V_{\text{Standard}}) \times (OD_{\text{Sample}(130\text{S})} - OD_{\text{Sample}(10\text{S})}) / (OD_{\text{Standard}} - \\ &\quad OD_{\text{Blank}}) / (V_{\text{Sample}} \times N / V_{\text{Assay}}) / T \\ &= 4000 \times (OD_{\text{Sample}(130\text{S})} - OD_{\text{Sample}(10\text{S})}) / (OD_{\text{Standard}} - OD_{\text{Blank}}) / N \end{aligned}$$

4. According to the volume of sample

$$\begin{aligned} \text{G6PDH (U/ml)} &= (C_{\text{Standard}} \times V_{\text{Standard}}) \times (OD_{\text{Sample}(130\text{S})} - OD_{\text{Sample}(10\text{S})}) / (OD_{\text{Standard}} - \\ &\quad OD_{\text{Blank}}) / V_{\text{Sample}} / T \\ &= 4000 \times (OD_{\text{Sample}(130\text{S})} - OD_{\text{Sample}(10\text{S})}) / (OD_{\text{Standard}} - OD_{\text{Blank}}) \end{aligned}$$

C_{Standard} : the standard concentration, 400 $\mu\text{mol/L}$ = 400 nmol/ml;

V_{Standard} : the volume of standard, 200 μl = 0.2 ml;

C_{Protein} : the protein concentration, mg/ml;

W: the weight of sample, g;

N: the quantity of cell or bacteria, $N \times 10^4$;


V_{Sample} : the volume of sample, 0.01 ml;

V_{Assay} : the volume of Assay buffer, 1 ml;


T: the reaction time, 2 minutes.

VII. TYPICAL DATA

The standard curve is for demonstration only. A standard curve must be run with each assay.

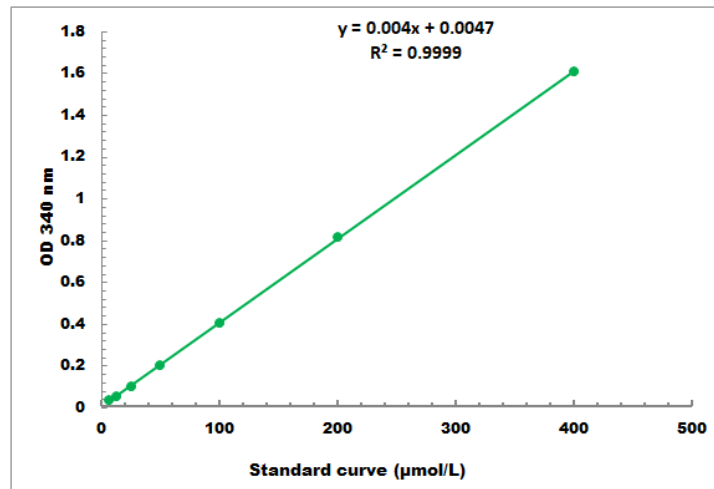

Detection Range: 4 µmol/L - 400 µmol/L

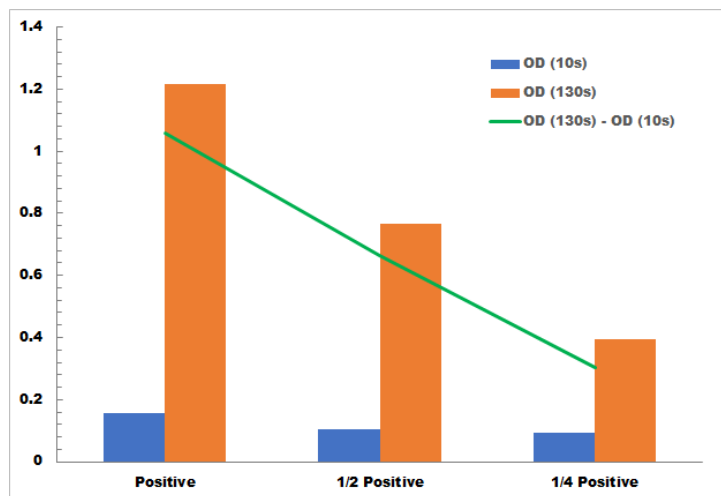

Positive Control reaction in 96-well plate assay with decreasing the concentration

VIII. TECHNICAL SUPPORT

For troubleshooting, information or assistance, please go online to www.cohesionbio.com or contact us at techsupport@cohesionbio.com

IX. NOTES