



# Beta-Glucuronidase Microplate Assay Kit User Manual

Catalog # FTA0151

(Version 1.2A)

Detection and Quantification of Beta-Glucuronidase (GGT) Activity in Urine, Serum, Plasma, Tissue extracts, Cell lysate, Cell culture media and Other biological fluids Samples.

For research use only. Not for diagnostic or therapeutic procedures.



I. INTRODUCTION	2
II. KIT COMPONENTS	3
III. MATERIALS REQUIRED BUT NOT PROVIDED	3
IV. SAMPLE PREPARATION	4
V. ASSAY PROCEDURE	5
VI. CALCULATION	6
VII. TYPICAL DATA	7
VIII. TECHNICAL SUPPORT	7
IX NOTES	7



### I. INTRODUCTION

Beta-Glucuronidase is hydrolytic enzyme responsible for the breakdown of carbohydrates. Specifically, beta-Glucuronidase cleave the terminal beta-D-glucuronic acid residue from the non-reducing terminus of a mucopolysaccharide chain. In humans, these enzymes are found in the lysosome of many tissue types. Loss of Beta-Glucuronidase activity results in a metabolic disease known as Sly syndrome. One pharmaceutical application for these enzymes is the metabolism of glucuronidated prodrugs into active pharmacological compounds. As expression and activities of beta-Glucuronidase vary substantially between tissue types and disease states, these enzymes have been used to achieve targeted activation of oncotherapeutic compounds, some of which may be toxic to healthy cells not associated with malignancy or disease. It is thus important to have knowledge of the beta-Glucuronidase activity in the tested sample to determine whether the prodrug or active form will predominate.

Beta-Glucuronidase Microplate Assay Kit is designed to measure mannitol in various samples in 96-well microplate. The color intensity at 560nm is directly proportional to beta-Glucuronidase activity in the sample.



### **II. KIT COMPONENTS**

Component	Volume	Storage	
96-Well Microplate	1 plate		
Assay Buffer	30 ml x 4	4 °C	
Reaction Buffer	10 ml x 1	4 °C	
Substrate	Powder x 1	4 °C	
Dye Reagent	10 ml x 1	4 °C	
Standard	Powder x 1	4 °C	
Plate Adhesive Strips	3 Strips		
Technical Manual	1 Manual		

# Note:

Substrate: add 2 ml Reaction Buffer to dissolve before use.

Standard : add 1 ml ethanol to dissolve before use, then add 20  $\mu l$  in to 980  $\mu l$ 

ethanol, the concentration will be 1000  $\mu$ mol/L.

# III. MATERIALS REQUIRED BUT NOT PROVIDED

- 1. Microplate reader to read absorbance at 560 nm
- 2. Distilled water
- 3. Pipettor, multi-channel pipettor
- 4. Pipette tips
- 5. Mortar
- 6. Ice
- 7. Centrifuge
- 8. Timer
- 9. Ethanol



# IV. SAMPLE PREPARATION

1. For tissue samples

Weigh out  $0.1\,\mathrm{g}$  tissue, homogenize with 1 ml Assay buffer on ice, centrifuged at  $10,000\mathrm{g}$  4 °C for 10 minutes, take the supernatant into a new centrifuge tube and keep it on ice for detection.

2. For serum, plasma and other biological fluids samples Detect directly.



### V. ASSAY PROCEDURE

Add following reagents in the microplate:

Reagent	Sample	Control	Standard	Blank		
Sample	10 μΙ					
Standard			10 μΙ			
Distilled water		10 μΙ	20 μΙ	30 μΙ		
Reaction Buffer	70 μl	70 μΙ	70 μl	70 μΙ		
Substrate	20 μΙ	20 μΙ				
Mix, then put it in the oven, incubate at 37 °C for 18 hours.						
Dye Reagent	100 μΙ	100 μΙ	100 μΙ	100 μΙ		
Mix, record absorbance measured at 560 nm.						

### Note:

- 1) Perform 2-fold serial dilutions of the top standards to make the standard curve.
- 2) For unknown samples, we recommend doing a pilot experiment & testing several doses to ensure the readings are within the standard curve range. If the enzyme activity is lower, please add more sample into the reaction system; or increase the reaction time; if the enzyme activity is higher, please dilute the sample, or decrease the reaction time.



### VI. CALCULATION

**Unit Definition:** One unit of Beta-Glucuronidase activity is the enzyme that produce 1 nmol of the Phenolphthalein acid per hour.

1. According to the protein concentration of sample

$$\begin{aligned} \text{Beta-GD (U/mg)} &= \left( \text{C}_{\text{Standard}} \times \text{V}_{\text{Standard}} \right) \times \left( \text{OD}_{\text{Sample}} - \text{OD}_{\text{Control}} \right) / \left( \text{OD}_{\text{Standard}} - \text{OD}_{\text{Blank}} \right) / \\ & \left( \text{V}_{\text{Sample}} \times \text{C}_{\text{Protein}} \right) / \text{T} \\ &= 55.56 \times \left( \text{OD}_{\text{Sample}} - \text{OD}_{\text{Control}} \right) / \left( \text{OD}_{\text{Standard}} - \text{OD}_{\text{Blank}} \right) / \text{C}_{\text{Protein}} \end{aligned}$$

2. According to the weight of sample

Beta-GD (U/g) = 
$$(C_{Standard} \times V_{Standard}) \times (OD_{Sample} - OD_{Control}) / (OD_{Standard} - OD_{Blank}) / (W$$

$$\times V_{Sample} / V_{Assay}) / T$$

$$= 55.56 \times (OD_{Sample} - OD_{Control}) / (OD_{Standard} - OD_{Blank}) / W$$

3. According to the volume of sample

Beta-GD (U/mI) = 
$$(C_{Standard} \times V_{Standard}) \times (OD_{Sample} - OD_{Control}) / (OD_{Standard} - OD_{Blank}) / V_{Sample} / T$$

$$= 55.56 \times (OD_{Sample} - OD_{Control}) / (OD_{Standard} - OD_{Blank})$$

 $C_{Standard}$ : the concentration of Standard, 1000 µmol/L = 1000 nmol/ml;

C<sub>Protein</sub>: the protein concentration, mg/ml;

W: the weight of sample, g;

V<sub>Sample</sub>: the volume of sample, 0.01 ml;

V<sub>Standard</sub>: the volume of standard, 0.01 ml;


V<sub>Assay</sub>: the volume of Assay buffer, 1 ml;

T: the reaction time, 18 hours.



# VII. TYPICAL DATA

The standard curve is for demonstration only. A standard curve must be run with each assay.

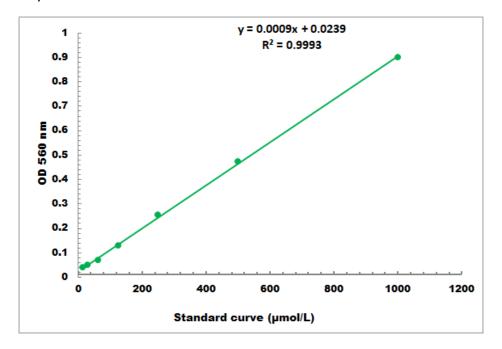

Detection Range: 10 μmol/L - 1000 μmol/L

### VIII. TECHNICAL SUPPORT

For troubleshooting, information or assistance, please go online to www.cohesionbio.com or contact us at techsupport@cohesionbio.com

# IX. NOTES