



# Histamine Microplate Assay Kit User Manual

Catalog # FTA0208

(Version 1.2A)

Detection and Quantification of Histamine Content in Tissue extracts and Other biological fluids Samples.

For research use only. Not for diagnostic or therapeutic procedures.



I. INTRODUCTION	2
II. KIT COMPONENTS	3
III. MATERIALS REQUIRED BUT NOT PROVIDED	3
IV. SAMPLE PREPARATION	4
V. ASSAY PROCEDURE	5
VI. CALCULATION	6
VII. TYPICAL DATA	7
VIII. TECHNICAL SUPPORT	7
IX NOTES	7



## I. INTRODUCTION

Histamine is a biogenic amine naturally present in many foods and body cells in amounts without toxicological significance. It is also a contaminant that mostly found in the Scombridae family fishes such as tuna and mackerel or other seafood products when improperly handled or stored. The consumption of foods containing high level of histamine may lead to an allergy-like food poisoning known as scombroid poisoning.

Histamine Microplate Assay Kit is a sensitive assay for determining Histamine content in various samples. The intensity of product color, measured at 480 nm is directly proportional to histamine concentration in the sample.



## **II. KIT COMPONENTS**

Component	Volume	Storage
96-Well Microplate	1 plate	
Assay Buffer	30 ml x 4	4 °C
Reaction Buffer	9 ml x 1	4 °C
Dye Reagent A	Powder x 1	4 °C
Dye Reagent B	Powder x 1	4 °C
Dye Reagent B Diluent	1 ml x 1	4 °C
Standard	Powder x 1	4 °C
Technical Manual	1 Manual	

## Note:

**Dye Reagent A:** add 8 ml Distilled water to dissolve before use.

Dye Reagent B: add 1 ml Dye Reagent B Diluent to dissolve before use.

**Dye Reagent Working Solution:** add all Dye Reagent B into Dye Reagent A, mix.

Standard : add 1 ml Assay Buffer to dissolve before use; then add 200  $\mu l$  into 800  $\mu l$ 

Assay Buffer, the concentration will be 2 mmol/L.

# III. MATERIALS REQUIRED BUT NOT PROVIDED

- 1. Microplate reader to read absorbance at 480 nm
- 2. Distilled water
- 3. Pipettor, multi-channel pipettor
- 4. Pipette tips
- 5. Mortar
- 6. Centrifuge
- 7. Timer



## IV. SAMPLE PREPARATION

# 1. For tissue samples

Weigh out 0.1 g tissue, homogenize with 1 ml Assay buffer on ice; incubate at 60 °C water bath for 30 minutes; centrifuged at 8000g 4 °C for 10 minutes, take the supernatant into a new centrifuge tube for detection.

# 3. For liquid samples

Detect directly, or dilute with Assay Buffer.



# V. ASSAY PROCEDURE

Warm all regents to room temperature before use.

Add following reagents into the microplate:

Reagent	Sample	Standard	Blank		
Reaction Buffer	90 μΙ	90 μΙ	90 μΙ		
Sample	20 μΙ				
Standard		20 μΙ			
Distilled water			20 μΙ		
Mix.					
Dye Reagent Working Solution	90 μΙ	90 μΙ	90 μΙ		
Mix, incubate at room temperature for 10 minutes, measured at 480 nm and record					
the absorbance.					

## Note:

- 1) Perform 2-fold serial dilutions of the top standards to make the standard curve.
- 2) The concentrations can vary over a wide range depending on the different samples. For unknown samples, we recommend doing a pilot experiment & testing several doses to ensure the readings are within the standard curve range.



# VI. CALCULATION

# 1. According to the volume of sample

Histamine (
$$\mu$$
mol/ml) = ( $C_{Standard} \times V_{Sample}$ ) × ( $OD_{Sample} - OD_{Blank}$ ) / ( $OD_{Standard} - OD_{Blank}$ )/
$$V_{Standard}$$
= 2 × ( $OD_{Sample} - OD_{Blank}$ ) / ( $OD_{Standard} - OD_{Blank}$ )

# 2. According to the weight of sample

Histamine (
$$\mu$$
mol/g) = ( $C_{Standard} \times V_{Sample}$ ) × ( $OD_{Sample} - OD_{Blank}$ ) / ( $OD_{Standard} - OD_{Blank}$ )

C<sub>Protein</sub>: the protein concentration, mg/ml;

 $C_{Standard}$ : the standard concentration, 2 mmol/L = 2  $\mu$ mol/ml;

W: the weight of sample, g;

V<sub>Sample</sub>: the volume of sample, 0.02 ml;


V<sub>Standard</sub>: the volume of standard, 0.02 ml;

V<sub>Assay</sub>: the volume of Assay buffer, 1 ml.



# VII. TYPICAL DATA

The standard curve is for demonstration only. A standard curve must be run with each assay.

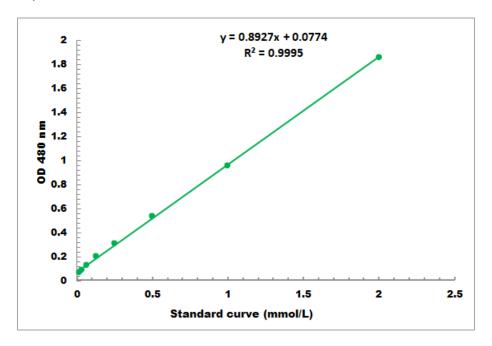

Detection Range: 0.02 mmol/L - 2 mmol/L

# VIII. TECHNICAL SUPPORT

For troubleshooting, information or assistance, please go online to www.cohesionbio.com or contact us at techsupport@cohesionbio.com

# IX. NOTES