

Goat Anti-PPP1R12A Antibody

Peptide-affinity purified goat antibody

Catalog # AF4330a


Product Information

Application	IF, FC, E
Other Accession	NP_002471.1 , NP_001231919.1 , NP_001231921.1 , NP_001137358.1
Reactivity	Human
Predicted	Human
Host	Goat
Clonality	Polyclonal
Isotype	IgG

Additional Information


Dilution	IF~~1:50~200 FC~~1:10~50 E~~N/A
Format	0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin.
Immunogen	Peptide with sequence C-REDEYKQKYSRTYD, from the C Terminus of the protein sequence according to NP_002471.1, NP_001231919.1, NP_001231921.1, NP_001137358.1.
Storage	Maintain refrigerated at 2-8°C for up to 6 months. For long term storage store at -20°C in small aliquots to prevent freeze-thaw cycles.
Precautions	Goat Anti-PPP1R12A Antibody is for research use only and not for use in diagnostic or therapeutic procedures.

Images

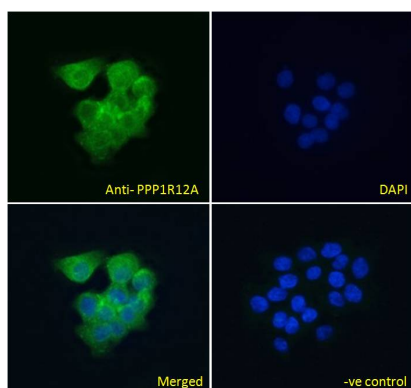


AF4330a Immunofluorescence analysis of paraformaldehyde fixed A431 cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (4ug/ml), showing cytoplasmic staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (4ug/ml).

AF4330a Immunofluorescence analysis of

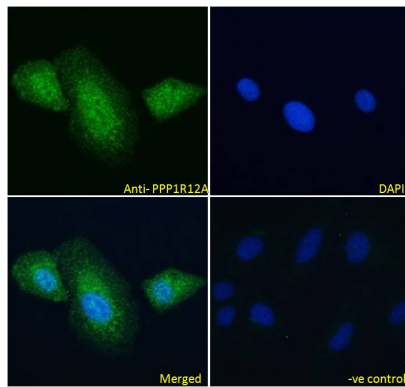

paraformaldehyde fixed U2OS cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (4ug/ml), showing cytoplasmic and some nuclear staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (4ug/ml)

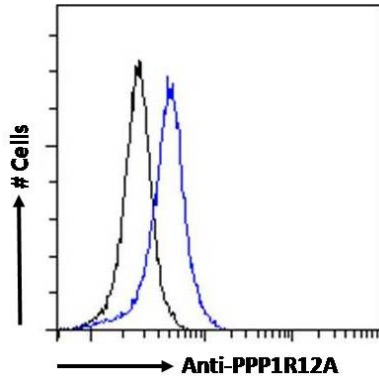

AF4330a Flow cytometric analysis of paraformaldehyde fixed HEK293 cells (blue line), permeabilized with 0.5% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (0.4ug/ml). IgG control: Unimmunized goat IgG (black line) followed by Alexa Fluor 488 secondary antibody.

Please note: All products are 'FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC OR THERAPEUTIC PROCEDURES'.