

Recombinant Anti-IKK gamma Rabbit mAb

Catalog # AP95041

Product Information

Application	WB, IHC, FC, IP, IF/IC
Primary Accession	Q9Y6K9
Reactivity	Human, Mouse, Rat
Host	Rabbit
Clonality	Monoclonal
Calculated MW	48198

Additional Information

Gene ID	8517
Other Names	FIP3; NEMO; NF-kappa-B essential modulator; NEMO; FIP-3; IκB kinase-associated protein 1; IKKAP1; Inhibitor of nuclear factor kappa-B kinase subunit gamma; I-kappa-B kinase subunit gamma; IKK-gamma; IKKG; IκB kinase subunit gamma; NF-kappa-B essential modifier
Target/Specificity	KLH-conjugated synthetic peptide encompassing a sequence within human IKK gamma protein. The exact sequence is proprietary.
Dilution	WB~~WB (1/500 - 1/1000) IHC~~1:100~500 FC~~1:10~50 IP~~N/A IF/IC~~N/A
Format	Liquid in PBS, pH 7.3, 50% glycerol, 0.05% BSA, and 0.05% Proclin300.
Storage	Store at -20 °C.Stable for 12 months from date of receipt

Protein Information

Name	IKBKG (HGNC:5961)
Synonyms	FIP3, NEMO
Function	Regulatory subunit of the IKK core complex which phosphorylates inhibitors of NF-kappa-B thus leading to the dissociation of the inhibitor/NF-kappa-B complex and ultimately the degradation of the inhibitor (PubMed: 14695475 , PubMed: 20724660 , PubMed: 21518757 , PubMed: 9751060). Its binding to scaffolding polyubiquitin plays a key role in IKK activation by multiple signaling receptor pathways (PubMed: 16547522 , PubMed: 18287044 , PubMed: 19033441 , PubMed: 19185524 , PubMed: 21606507 , PubMed: 27777308 , PubMed: 33567255). Can recognize and bind both 'Lys-63'-linked and linear polyubiquitin upon cell stimulation, with a much higher affinity for linear polyubiquitin (PubMed: 16547522 , PubMed: 18287044 , PubMed: 19033441 , PubMed: 19185524 , PubMed: 21606507 , PubMed: 27777308). Could be implicated in NF-kappa-B-mediated protection

from cytokine toxicity. Essential for viral activation of IRF3 (PubMed:[19854139](#)). Involved in TLR3- and IFIH1-mediated antiviral innate response; this function requires 'Lys- 27'-linked polyubiquitination (PubMed:[20724660](#)).


Cellular Location

Cytoplasm. Nucleus Note=Sumoylated NEMO accumulates in the nucleus in response to genotoxic stress.


Tissue Location

Heart, brain, placenta, lung, liver, skeletal muscle, kidney and pancreas


Images

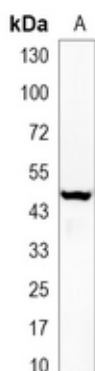

Western blot analysis of IKK gamma expression in Jurkat (A) whole cell lysates. (Predicted band size: 48 kD; Observed band size: 48 kD)

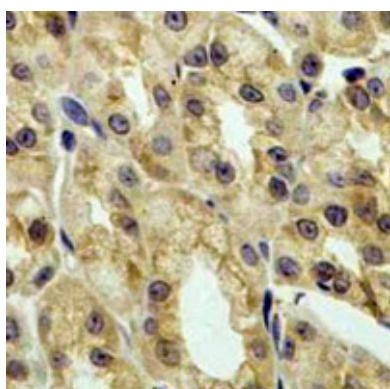

Immunohistochemical analysis of IKK gamma staining in mouse colon formalin fixed paraffin embedded tissue section. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH 6.0). The section was then incubated with the antibody at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

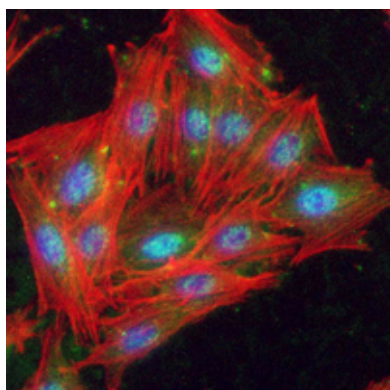

Immunofluorescent analysis of IKK gamma staining in Jurkat cells. Formalin-fixed cells were permeabilized with 0.1% Triton X-100 in TBS for 5-10 minutes and blocked with 3% BSA-PBS for 30 minutes at room temperature. Cells were probed with the primary antibody in 3% BSA-PBS and incubated overnight at 4 °C in a humidified chamber. Cells were washed with PBST and incubated with an AF488-conjugated secondary antibody (green) in PBS at room temperature in the dark. Phalloidin - AF594 was used to stain Actin filaments (red). DAPI was used to stain the cell nuclei (blue).

Please note: All products are 'FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC OR THERAPEUTIC PROCEDURES'.