

Anti-MOUSE IgG (H&L) (Peroxidase Conjugated) Secondary Antibody

Rabbit Polyclonal, Peroxidase (Horseradish)

Catalog # ASR1935

Product Information

Description	Anti-MOUSE IgG (H&L) (RABBIT) Antibody Peroxidase Conjugated
Host	Rabbit
Conjugate	Peroxidase (Horseradish)
Target Species	Mouse
Clonality	Polyclonal
Application	WB
Physical State	Lyophilized
Host Isotype	IgG
Target Isotype	IgG (H&L)
Buffer	0.02 M Potassium Phosphate, 0.15 M Sodium Chloride, pH 7.2
Immunogen	Mouse IgG whole molecule
Reconstitution Volume	1.0 mL
Reconstitution Buffer	Restore with deionized water (or equivalent)
Stabilizer	10 mg/mL Bovine Serum Albumin (BSA) - Immunoglobulin and Protease free
Preservative	0.01% (w/v) Gentamicin Sulfate. Do NOT add Sodium Azide!


Additional Information

Shipping Condition	Ambient
Application Note	When choosing a secondary antibody product, consideration must be given to species and immunoglobulin specificity, conjugate type, fragment and chain specificity, level of cross-reactivity, and host-species source and fragment
Purity	HRP Secondary Antibody Conjugate was prepared from monospecific antiserum by immunoaffinity chromatography using Mouse IgG coupled to agarose beads. Assay by immunoelectrophoresis resulted in a single precipitin arc against anti-Peroxidase, anti-Rabbit Serum, Mouse IgG and Mouse Serum.
Storage Condition	Store mouse second antibody conjugate at 4° C prior to restoration. For extended storage aliquot antibody and freeze at -20° C or below. Avoid cycles of freezing and thawing. Centrifuge product if not completely clear after standing at room temperature. This product is stable for several weeks at 4° C as an undiluted liquid. Dilute only prior to immediate use.
Precautions Note	This product is for research use only and is not intended for therapeutic or diagnostic applications.

Background

Mouse secondary antibody conjugated to horseradish peroxidase is available in a variety of formats. Anti-IgG Secondary Antibody conjugate is suitable for western blot, ELISA, immunoperoxidase electron microscopy and immunohistochemistry as well as other peroxidase-antibody based enzymatic assays requiring lot-to-lot consistency.

Images

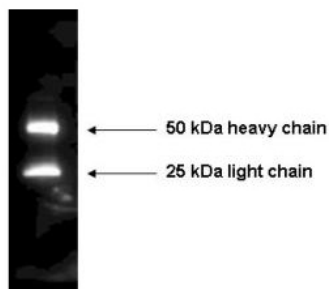

Western Blot of Abcepta's Anti-Mouse IgG Antibody Peroxidase Conjugated detecting mouse IgG. Mouse IgG (1.0 mg) was separated on a 4-20% gradient gel and then transferred to Whatman Protran BA 85 nitrocellulose (0.45 μ m pore size), and blocked with 5% BLOTTO (p/n B501-0500) o/n at 4° C. The membrane was then probed with a 1:2,500 dilution of HRP Rb-a-Mouse IgG [H&L] (p/n ASR1935) for 45 min at RT. FemtoMax Super Sensitive Chemiluminescent HRP Substrate (p/n FEMTOMAX-110). Was added and the image was captured using a BioSpectrum Imaging System with a BioChemi 500 -28C cooled 4.0 Megapixel CCD Camera from UVP (www.uvp.com, Upland, CA). Binning: 5X5. Exposure time: 2 min. Aperture open. No emission or excitation filters used.

Please note: All products are 'FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC OR THERAPEUTIC PROCEDURES'.